

UR free to go

Your portable urine collector

KANG CHIEN MEDICAL CO., LTD.

Where it comes from

It is very important for us to stay hydrate

70% of our body is water

Body's metabolism

Allowing free circulation



The IDEA

Harmful effects of dehydration in our body :

1. Painful kidney, ureter and urethra stones
2. Urinary bladder inflammation
3. Effect Heart rate, heart failure and kidney failure

The damaging effects of holding urine :

1. Frequent urination, especially in the middle of the night
2. Easily lead to urine stones or inflammation of the urinary system
3. Back Flow of the urine from bladder into kidney may destroy its function

Reason for URFTG

What elderly thinks:

1. Go for a bit of walk
2. Spend time with family and friends
3. Learn new things
4. JUST GO OUT



Reasons

However, in reality...



Others...

Things we see on the market these days.



Others...

And when we thought we could finally cover up the problem...

- *Diaper rash
- *Over flow
- *The Smells



New Product: UR free to go



URFTG :

- 1.Can wear daily
- 2.Easily clean
- 3.Reusable
- 4.Odorless

You are always free to go

Product



You can use a 20 CC pump to adjust for your own comfort.

- Air adjustable Valve:



Product



- Air vent :



Product



Users can use the same 20cc pump to injects water to the flushing hole for cleaning.

- Flushing hole :



Product



1. Soft material for long term wearing

2. Can adjust to fit the user's penis

- Air cushion :




Product




Soft inner ring for the user's penis

- Inner ring :

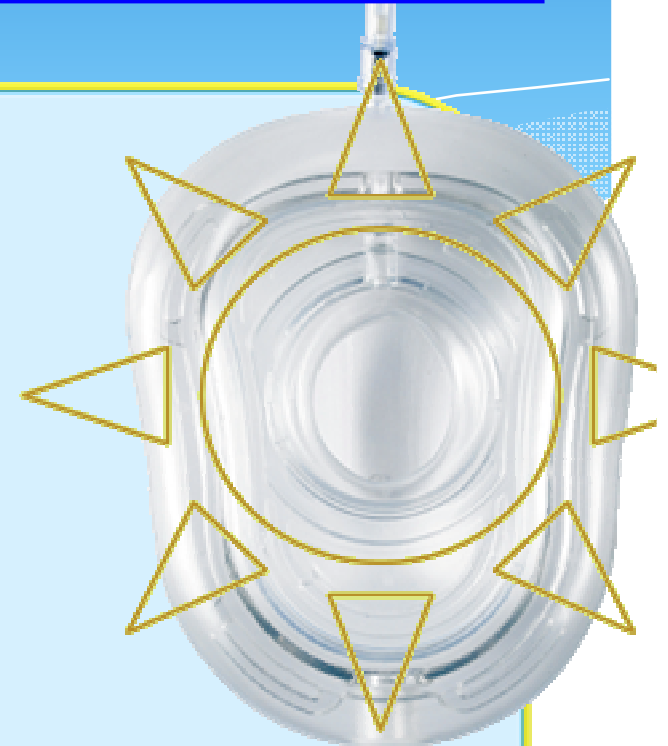

Product

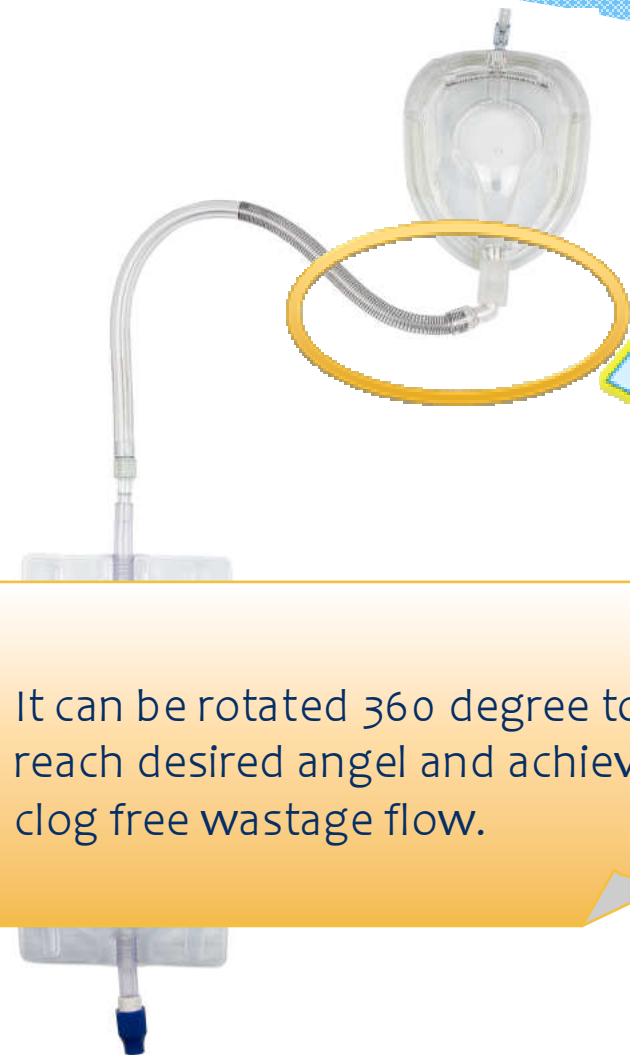

It can be rotated 360 degree to reach desired angel and achieve clog free wastage flow.

- **Metal wired tube with L connector**




Product



Connect to URFTG cup. Hold 500 cc

- Urine collection bag :

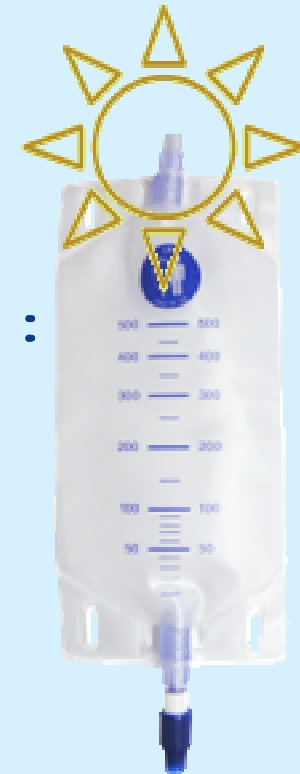

Product



Just twist the nob and let the unwanted waste go~



Product

- **Urine release nob:**



Product User Manual




How to wear the device



- Put the cup into the fixing underwear.



- Pull the tube out through the lower hole.

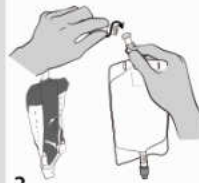

- Screw the tube firmly onto the Drip bag.



- Place the waist band around the waist and pull it firm.



- Pull the hip strap around and down the crotch, and secure it on both sides.



- Pull the inflation hole through the upper hole of the fixing underwear.




- Use the flushing syringe to inflate or deflate, to allow for appropriate space in relation to the individual genitals.



- Pull the strap through the two holes at the top of the Drip bag, with the anti-slip side inside, touching the skin. Adjust to a tightness that is comfortable and fixates the bag.

* Be careful not to damage the drip bag by too much pressure.

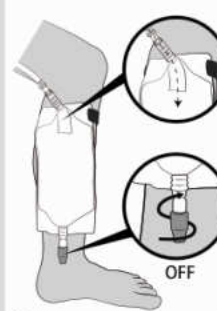

- Before use, make sure drip bag's outlet is closed. During use, make sure the inlet valve is open.



Complete chart.

Cleaning method

Easy flush while wearing




Flush the cup with clean water of less than 40°C, inserting it through the flushing inlet. During flushing the water will flow down the cup, it will not touch the body.

Thorough cleaning after taking it off

- From the support ring inject a large amount of water to wash the mask shape cup.
- Flush the water through the tube to the liquid collection bag. Then shake and squeeze the bag, to make the water thoroughly rinse every corner. Next, open the outlet valve to drain the water used to clean the device.


- You may use a soft dry cloth to gently wipe the inside of the cup. Then store the device where the temperature is between -20° and 45° C (-4° to 113° F).

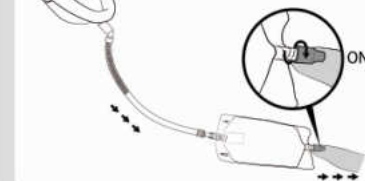


Product User Manual




Contents and the name of each part

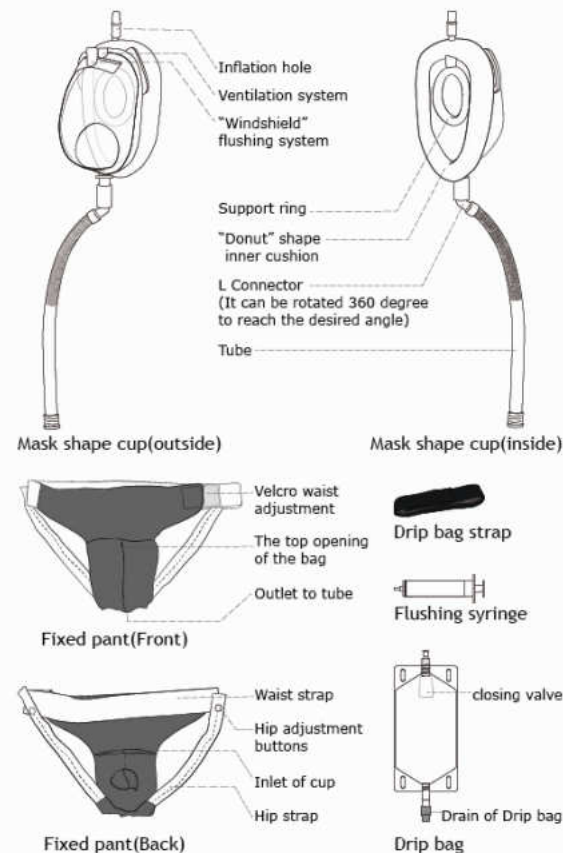

Points of attention

- 


This product can be cleaned and reused. A replacement each month is suggested.
- 

This product is designed to assist men who cannot go to the toilet in time, for example: long-haul truck drivers, those bothered by frequent urination or urinary incontinence, etc.
- 

This product is meant for individual use only. Do not share, to avoid transmitting diseases through contact.
- 

This product should be kept in a dry environment, and be hung out to dry after cleaning.
- 

The device can be washed using clean water of less than 40°C and agents suitable for washing the body. Do not use detergents, alcohol, or other corrosive agents.
- 

Store the device where the temperature is between -20° and 45° C (-4° to 113° F). Don't leave the device in your car, because temperatures in parked cars can exceed this range.
- 

When using this device in sitting position, put your leg and knee straight, to ensure smooth downward flow of the urine.
- 

This product is not suitable for use in lying position.

UR free to go Portable Urinal for Men User manual



Package Contents

Mask shape cup x 1
Drip bag x 2
Drip bag strap x 1
Fixed pant (optional)
Flushing syringe x 1
User manual x 1

Optional accessories with UR free to go

Fixed pant
Leg bag holder
Charcoal pant (Briefs/Boxers)
Drip bag replacements
(500/800/2000ml)

Product Lines

Product Name	Specs
<p>UR Free to go Portable Urinal for Men</p> 	<p>Mask-Shape Cup x1 Fixed Pantx1(optional) 500cc Drip bagx2 Drip bag strapx1 Flushing syringex1 User Manualx1</p>
<p>Fixed Pant</p> 	<p>L (28" ~ 42") XL (38" ~ 54")</p>
<p>Leg Bag Holder</p> 	<p>Free Size (28" ~ 42")</p>
<p>Charcoal Pant (Briefs/Boxers)</p> 	<p>L Size (28" ~ 42") XL (40" ~ 52) XXL(50" -62")</p>
<p>Drip Bag</p> 	<p>500 cc</p>
	<p>800 cc</p>
	<p>2,000 cc</p>

Target customers


1. Patient with urine problem
2. Patient with nerve damage , Anyone with mobility problems. Long term diaper users who want to improve the quality of Life
3. Radiation, dust-free working environment who must wear protective isolation contained work clothes.
4. Small aircraft, helicopters or gliders flying
5. Replace the urinal at night, avoid dark falling injuries.



Diapers VS Free to Go

Diaper	Free to GO
Need place to Change	Worn easily
6 piece of diaper each day	Washable. NO need to change up to 10 days
Need 2 people to change and clean.	Does not require help from others~ and cleaning only use 1 minute.
Takes 15 mines to clean and Change	Takes less than 10 minutes change to a new unit
Easy to cause Diaper rash and skin	NO diaper rash and easy on the skin
Terrible urine odor	NO odor
Leaks and cannot with take a large amount at one time	NO leaks and can with take large volume of Urine
Very Bulky and awkward looking	Less bulky and easy to conceal with in the users clothing.

Patent

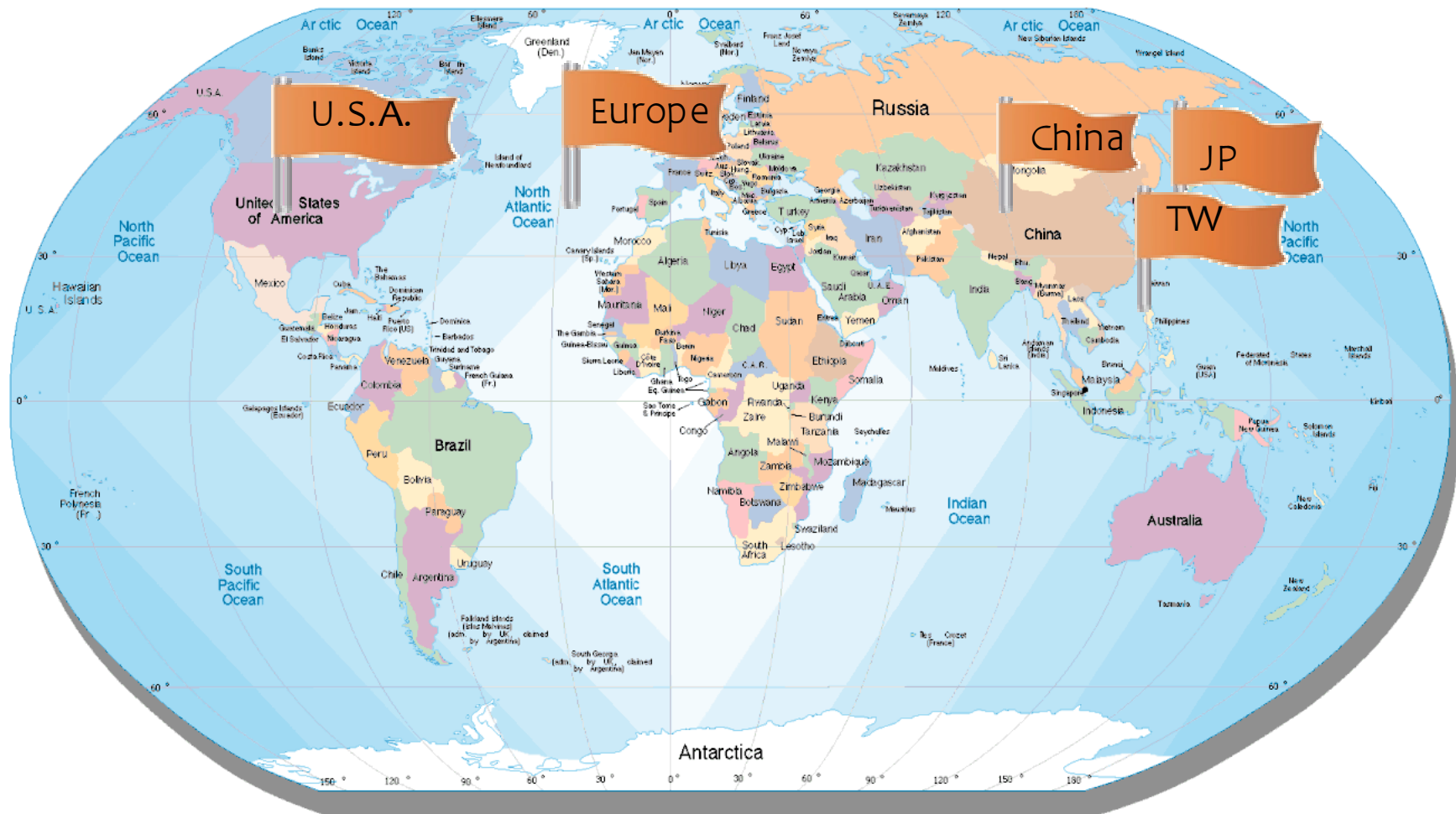

UR free to go Economic Benefits

Number of Days	US\$1.00 per pc and 6pcs a day	Kimberly Clark 《Depend》 Maximum Protection Brief with Easy Grip™ Tabs	UR Free to Go US\$99.99 (MSRP for one package)	Money You Save
1	6pcs	US\$6	US\$99.99	
30	180pcs	US\$180	US\$99.99	US\$80.01
60	360pcs	US\$360	US\$199.98 (US\$99.99 for another package)	US\$160.02
90	540pcs	US\$540	US\$299.97 (US\$99.99 for another package)	US\$240.03
180	1,080pcs	US\$1,080	US\$599.94 (US\$299.97 for another 3 packages)	US\$480.06
360	2,160pcs	US\$2,160	US\$1199.88 (US\$599.94 for another 6 packages)	US\$960.12 (45% of Saving)

SGS Report for Skin Irritation and Cytotoxicity Test

SGS

TEST REPORT

Mechanical & Hardgoods Lab

Report No. : HL80121/2011

Page : 1 of 4

Date : AUG. 16, 2011

Kang Chien Co. Ltd
9F.-2, No. 207, Dunhua N. Rd., Songshan Dist., Taipei 105, Taiwan

The following merchandise was submitted and identified by the applicant as:

Product Description: Free to go
Style/Item No.: KC00001
Manufacturer/Vendor: Foremount Enterprise CO. LTD
Country of Origin: Taiwan

We have tested the submitted sample(s) as requested and the following results were obtained:

Test Requested:
ISO 10993-5:2009 Cytotoxicity Test

Test Methods & Results:
--See following sheet(s)--

Date of Receipt: Aug. 05, 2011

Testing Period: Aug. 06, 2011 ~ Aug. 15, 2011

Conclusion: The submitted sample which the cell number for the test article extract was not over 50% less than the negative control and cell reactivity grading was equal to 0 (none). Therefore, the test article showed "negative" cytotoxicity.

Signed for and on behalf of
SGS Taiwan Ltd.


Benson Liao
Asst. Manager

Unless otherwise stated the results shown in this test report refer only to the sample(s) tested. This test report cannot be reproduced, except in full, without prior written permission of the Company. 除非另有說明，此報告內容僅針對測試樣品負責，未經本公司書面許可，不可部分複製。

This document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at www.sgs.com/taichung_and_conditions.htm and for electronic format documents, subject to Terms and Conditions for Electronic Documents at www.sgs.com/taichung_and_conditions.htm. Attention is drawn to the limitation of liability, indemnification and jurisdiction clauses defined therein. Any holder of this document is advised that information contained herein reflects the Company's findings at the time of its issuance only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not constitute parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law.

SGS Taiwan Ltd. | 125, 127, Wu Kang Road, Wu Ku District, New Taipei City, Taiwan / 新北市五股區五股工業區五工路125號, 127號
台灣檢驗科技股份有限公司 | 1886-21 2799-3279 | 1886-21 2299-2920 | www.tai.sgs.com

Member of SGS Group

SGS

TEST REPORT

Mechanical & Hardgoods Lab

Report No. : HL80122/2011

Page : 1 of 6

Date : AUG. 29, 2011

Kang Chien Co. Ltd
9F.-2, No. 207, Dunhua N. Rd., Songshan Dist., Taipei 105, Taiwan

The following merchandise was submitted and identified by the applicant as:

Product Description: Free to go
Style/Item No.: KC00001
Manufacturer/Vendor: Foremount Enterprise CO. LTD
Country of Origin: Taiwan

We have tested the submitted sample(s) as requested and the following results were obtained:

Test Requested:
ISO 10993-10:2010 Skin Irritation Test

Test Methods & Results:
--- see attached sheets---

Date of Receipt: Aug. 05, 2011

Testing Period: Aug. 08 ~ 26, 2011

Conclusion: No erythema and no oedema were observed on the skin of the rabbits. The Primary Irritation Index for the test article was calculated to be 0, and the overall response fell into the category "negligible", therefore the test article under this study showed "negative reaction".

Signed for and on behalf of
SGS Taiwan Ltd.


Benson Liao
Asst. Manager

Unless otherwise stated the results shown in this test report refer only to the sample(s) tested. This test report cannot be reproduced, except in full, without prior written permission of the Company. 除非另有說明，此報告內容僅針對測試樣品負責，未經本公司書面許可，不可部分複製。

This document is issued by the Company subject to its General Conditions of Service printed overleaf, available on request or accessible at www.sgs.com/taichung_and_conditions.htm and for electronic format documents, subject to Terms and Conditions for Electronic Documents at www.sgs.com/taichung_and_conditions.htm. Attention is drawn to the limitation of liability, indemnification and jurisdiction clauses defined therein. Any holder of this document is advised that information contained herein reflects the Company's findings at the time of its issuance only and within the limits of Client's instructions, if any. The Company's sole responsibility is to its Client and this document does not constitute parties to a transaction from exercising all their rights and obligations under the transaction documents. This document cannot be reproduced except in full, without prior written approval of the Company. Any unauthorized alteration, forgery or falsification of the content or appearance of this document is unlawful and offenders may be prosecuted to the fullest extent of the law.

SGS Taiwan Ltd. | 125, 127, Wu Kang Road, Wu Ku District, New Taipei City, Taiwan / 新北市五股區五股工業區五工路125號, 127號
台灣檢驗科技股份有限公司 | 1886-21 2299-3279 | 1886-21 2299-2920 | www.tai.sgs.com

Member of SGS Group

USA FDA Registered Number

- * D-U-N-S Number: 656944027
- * Business Name: Kang chien Medical Co Ltd
- * City: Taipei City
- * Country: Taiwan
- * 9f.-2, No.207, Dunhua N. Rd., Songshan Dist., Taipei City 10595, Taiwan (R.O.C.)
- * 876.5250..... Urine collector and accessories
- * Owner Operator Name: 10045755

USA Medicare Coverage Number

“UR FREE TO GO” Portable Urinal for Men (Small Cup)	Part Number KC-001A	HCPCS Code A5105
“UR FREE TO GO” Portable Urinal for Men (Big Cup)	Part Number KC-0002	HCPCS Code A5105
“UR FREE TO GO” Drip Bag Replacement (500ML)	Part Number KCU-01A5	HCPCS Code A4358
“UR FREE TO GO” Drip Bag Replacement (800ML)	Part Number KCU-01A8	HCPCS Code A4358
“UR FREE TO GO” Drip Bag Replacement (2000ML)	Part Number KCU-01A20	HCPCS Code A4358

Product insurance (liability)



Fubon Insurance Co. Ltd.
237, Chen Kuo South Road, Sec. 1, Taipei, Taiwan, Republic of China
TEL: 886-2-2706 7890
免費申訴電話 (Toll-Free Number) : 0800-009-888
100.09.15 富保業字第 1000001429 號函審查

Commercial General Liability Insurance (Claim-Made Form)
PRODUCTS LIABILITY COVERAGE
POLICY NO.: 0500-02PD000039

ORIGINAL

Declaration

Named Insured & Mailing Address:

Kang Chien Co., Ltd.
9F -2, No.207, Dunhua N. Rd., Songshan Dist., Taipei City 105, Taiwan (R.O.C.)

Applicant: Kang Chien Co., Ltd.

Policy Period: From February 01, 2013 to February 01, 2014 (from noon of standard time at your mailing address shown above)

Retroactive Date: February 01, 2013

Extended Reporting Periods: Nil

Business Description: Manufacturer

Insured Products: Portable Urinal Collector (Excluding Products designed for Infant), which are manufactured / distributed by the insured and sold to the Coverage Territory of this policy.
Only the products sold via the vendors on Vendors List are insured under this policy.
(Warranted that all products insured are complied with statutory requirement and other mandatory and / or voluntary industrial / safety standards of importing countries.)

Insured Premises: Not Applicable

Limits of Insurance: (Supplementary Payment will be included within Limits of insurance)

General Aggregate Limit	USD 1,000,000.-
Coverage- A Premises – Operation Liability Aggregate Limit	Not Covered
Products Liability Aggregate Limit	USD 1,000,000.-
Completed Operations Aggregate Limit	Not Covered
Coverage- B Personal & Advertising Injury Aggregate Limit	Not Covered
Each Occurrence Limit	USD 1,000,000.-
Damage to Premises Rented to You Limit	Not Covered

Deductible: (Legal cost and expenses inclusive)
USD 2,500.- each occurrence

Estimated Annual Turnover: USA / Canada : USD 67,000.-
Rest of World : USD 100,000.- (Est. Turnover including Taiwan)

Adjustable Rate: 2.4 o/o

Minimum / Deposit Premium: USD 3,500.-
(Actual annual premium will be adjusted based on the actual sales at the end of policy year);

Coverage Territory / Jurisdiction: Worldwide including USA / Canada / Taiwan
(However, under no circumstances shall any country or jurisdiction which is the subject of any sanction, prohibition or restriction under United Nations resolutions or the trade or economic sanctions, laws or regulations of the European Union, United Kingdom or United States of America is being covered under this policy.)

The authority of the Insurer has hereunto set his hand at Taipei, this 21st day of Jan., 2013.

經理人 陳燦煌

AWARDS

GENEVA Golden award



AWARDS

Malaysia ITEX award



AWARDS

Taiwan invention award



Wish

When you need help on walking,
you use...

SO~ Why not?



Expectations

One small change of habit will result in a great improvement of life!



Thank you

North America Sale and After Service Office

1624 West 75th Avenue, Vancouver,
BC, V6P 6G2 Canada

Direct Line: 778-800-9406

Fax Number: 604-263-9529

Email: john.urftg@gmail.com

Website: www.urftg.com